


For short implant placement, implant major diameter needs to be $\leq$ the bur (average diameter) at the 8mm laser mark. Please refer to page 16 in the Instructions for Use Manual.

In Ridge Expansion cases, please oversize your osteotomy and make sure that the crest diameter is equal to or larger than the implant major diameter.

Use Densah Burs in full-step increments for Sinus Lift cases. Example: 2.0mm, 3.0mm, 4.0mm, 5.0mm

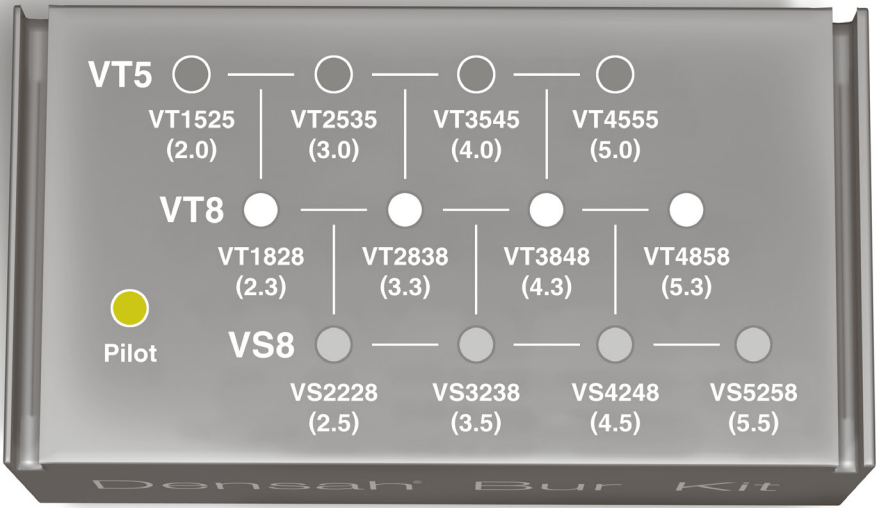

Use large block display to compare Bur identification system when using the schematic below for proper Bur usage

● VT5 Set ○ VT8 Set ● VS8 Set

Densifying Mode CCW (800-1500) RPMs / Cutting Mode CW (800-1500) RPMs																			
Dentium			SuperLine, Implantium, and SuperLine II																
			Soft Bone							Hard Bone (Mandible)									
										In densifying mode make sure your osteotomy is 1.0 mm deeper than the actual implant final length. In extreme hard bone, utilize DAC (Densify After Cut) Protocol. Find protocol in IFU.									
Geometry	Major Ø	Minor Ø	Pilot	Bur 1	Bur 2	Bur 3	Bur 4	Bur 5	Densah® Bur Block Display	Pilot	Bur 1	Bur 2	Bur 3	Bur 4	Bur 5	Bur 6	Bur 7	Bur 8	Densah® Bur Block Display
Slow Taper	3.4		Pilot	VT1828 (2.3)	VS2228* (2.5)	—	—			Pilot	VT1525 (2.0)	VT1828 (2.3)	VT2535* (3.0)	—	—	—	—	—	
Slow Taper	3.8		Pilot	VT1525 (2.0)	VT2535* (3.0)	—	—			Pilot	VT1828 (2.3)	VT2535 (3.0)	VT2838* (3.3)	—	—	—	—	—	
Slow Taper	4.3		Pilot	VT1828 (2.3)	VT2838 (3.3)	VS3238* (3.5)	—			Pilot	VT1828 (2.3)	VT2535 (3.0)	VT2838 (3.3)	VT3545* (4.0)	—	—	—	—	
Slow Taper	4.8		Pilot	VT1525 (2.0)	VT2535 (3.0)	VT3545* (4.0)	—			Pilot	VT1525 (2.0)	VT2535 (3.0)	VT2838 (3.3)	VT3545 (4.0)	VT3848* (4.3)	—	—	—	
Straight	4.8		Pilot	VT1828 (2.3)	VT2838 (3.3)	VT3848* (4.3)	—			Pilot	VT1828 (2.3)	VT2535 (3.0)	VT2838 (3.3)	VT3545 (4.0)	VT3848 (4.3)	VS4248* (4.5)	—	—	
Slow Taper	5.8		Pilot	VT1525 (2.0)	VT2535 (3.0)	VT3545 (4.0)	VT4555* (5.0)			Pilot	VT1828 (2.3)	VT2535 (3.0)	VT2838 (3.3)	VT3545 (4.0)	VT3848 (4.3)	VT4858* (5.3)	—	—	
Taper	6.0	4.8	Pilot	VT1828 (2.3)	VT2838 (3.3)	VT3848 (4.3)	VT4858* (5.3)			Pilot	VT1828 (2.3)	VT2535 (3.0)	VT2838 (3.3)	VT3545 (4.0)	VT3848 (4.3)	VT4858* (5.3)	VS5258* (5.5)	—	
Taper	7.0	5.8	Pilot	VT1828 (2.3)	VT2838 (3.3)	VT3848 (4.3)	VT4858 (5.3)	VS5258* (5.5)		Pilot	VT1828 (2.3)	VT2535 (3.0)	VT2838 (3.3)	VT3545 (4.0)	VT3848 (4.3)	VT4858 (5.3)	VS5258* (5.5)	—	


*Denotes implant placement.

*Clinician experience and judgment should be used in conjunction with the Densifying Reference Guide recommendation and suggested use protocols.

For short implant placement, implant major diameter needs to be $\leq$ the bur (average diameter) at the 8mm laser mark. Please refer to page 16 in the Instructions for Use Manual.

In Ridge Expansion cases, please oversize your osteotomy and make sure that the crest diameter is equal to or larger than the implant major diameter.

Use Densah Burs in full-step increments for Sinus Lift cases. Example: 2.0mm, 3.0mm, 4.0mm, 5.0mm

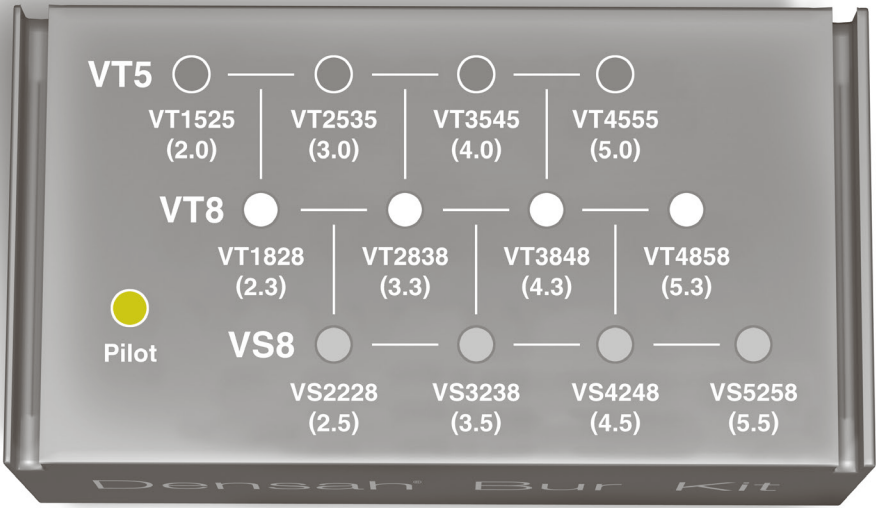

Use large block display to compare Bur identification system when using the schematic below for proper Bur usage

● VT5 Set ○ VT8 Set ● VS8 Set

Densifying Mode CCW (800-1500) RPMs / Cutting Mode CW (800-1500) RPMs

Dentium			Implant II														
			Soft Bone						Hard Bone (Mandible)								
									In densifying mode make sure your osteotomy is 1.0 mm deeper than the actual implant final length. In extreme hard bone, utilize DAC (Densify After Cut) Protocol. Find protocol in IFU.								
Geometry	Major Ø	Minor Ø	Pilot	Bur I	Bur 2	Bur 3	Bur 4	Densah® Bur Block Display	Pilot	Bur I	Bur 2	Bur 3	Bur 4	Bur 5	Bur 6	Bur 7	Densah® Bur Block Display
Slow Taper	3.6		Pilot	VTI525 (2.0)	VS2228* (2.5)	—	—		Pilot	VTI525 (2.0)	VTI828 (2.3)	VT2535* (3.0)	—	—	—	—	
Slow Taper	4.0		Pilot	VTI828 (2.3)	VT2838* (3.3)	—	—		Pilot	VTI828 (2.3)	VT2535 (3.0)	VT2838 (3.3)	VS3238* (3.5)	—	—	—	
Slow Taper	4.5		Pilot	VTI525 (2.0)	VT2535 (3.0)	VT3545* (4.0)	—		Pilot	VTI525 (2.0)	VT2535 (3.0)	VT2838 (3.3)	VT3545 (4.0)	VT3848* (4.3)	—	—	
Slow Taper	5.0		Pilot	VTI828 (2.3)	VT2838 (3.3)	VT3848* (4.3)	—		Pilot	VTI828 (2.3)	VT2535 (3.0)	VT2838 (3.3)	VT3545 (4.0)	VT3848 (4.3)	VS4248* (4.5)	—	

*Denotes implant placement.

*Clinician experience and judgment should be used in conjunction with the Densifying Reference Guide recommendation and suggested use protocols.